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

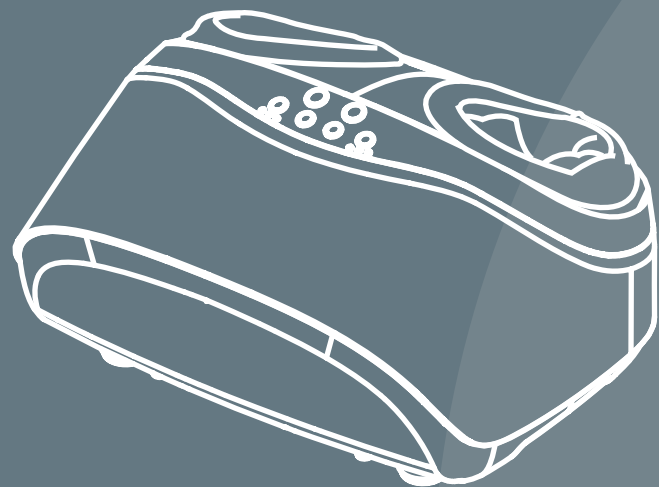
## IMPORTANT SAFETY INSTRUCTIONS

### OPERATING SURROUNDINGS

- Do not use or store the foot massager under damp environment: i.e. bathroom, sauna, spa to avoid any moldy, electric shock, fire or mechanical failure.
- Don't use the foot massager outdoors
- Don't use the foot massager under dusty, smeary and anoxic conditions
- Don't place the foot massager near heater, stove or sunlight radiation to avoid inner components overheating to cause fire
- Make sure to place the foot massager on a flat level ground
- The bracket part should be placed on the plain ground
- To avoid an electric shock, please make sure to connect an earth wire
- Don't use the machine while smoking
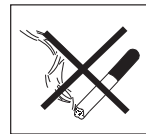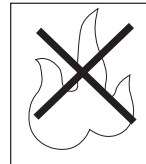
### FOLLOWING CONDITIONS SHOULD CONSULT DOCTORS BEFORE USE

- Suffered from malignant tumor
- Serious heart disease
- Serious skin disease
- Serious boned hyperplasia
- In pregnancy or during the menses
- Be on the rest cure under the doctor's supervision
- The back bended exceptionally by the accident or disease
- Needing peaceful surroundings
- Arthritis disease
- Implanted pacemaker for medical purpose
- Suffered from bowels disease (i.e. gastritis, enteritis, diarrhea, hepatitis)
- Serious blood circulation trouble
- Consciousness trouble
- Be in hyper pyretic

# IMPORTANT SAFETY INSTRUCTIONS

For safety reasons and in order to clearly indicate potential dangers should the product be used improperly, instructions have been marked in the following manners and should be strictly noted

## ⊘ CAUTIONS

**Non-professional technicians are not allowed to dismantle or repair to avoid fire, electrical shock or serious hurt**

**Do not flush to the switch when cleaning to avoid insulation fault, electrical shock or short-circuit**

**When storing or maintenance, be sure to disconnect the power plug from socket**

Do not detach the plug by wet hands

**Make sure to shut off and detach the power plug when cleaning or storing product**

During cleaning or moving the machine, make sure to detach power cord to avoid fire or electrical shock

**Do not place the machine onto an unstable surface to avoid roast, insulation failure, electrical shock or turnover**

When storing the machine, make sure to disconnect power plug to avoid fire caused by electricity leaking

**Clean power plug periodically to avoid fire due to dust**

**Do not use it at any unsafe places or on other purposes**

It is the main reason of product damage or accident

**Do not store the machine onto an unstable surface**

Unexpected earthquake or external force will cause damage or injury




**Keep supply cord unobstructed at all time, i.e. not twisted, knotted, over flexed or pulled and do not place anything heavy on it**

# PARTS NAME


NAME OF EACH PART

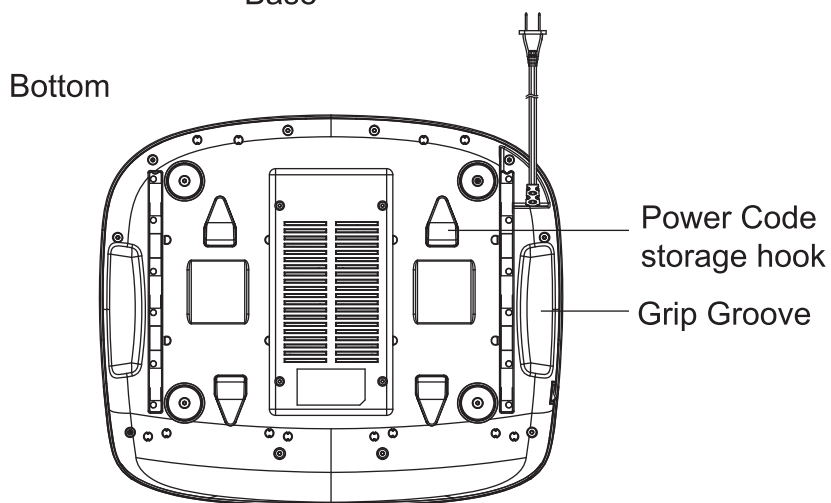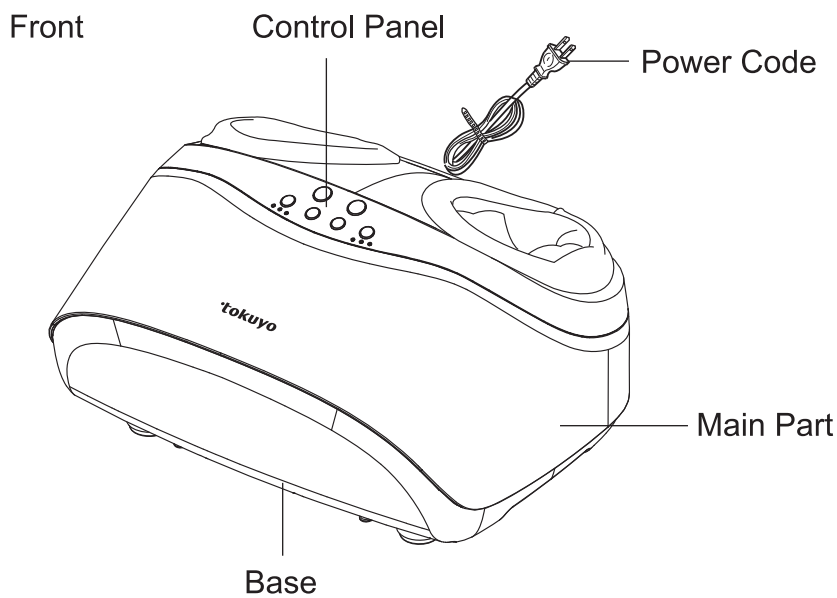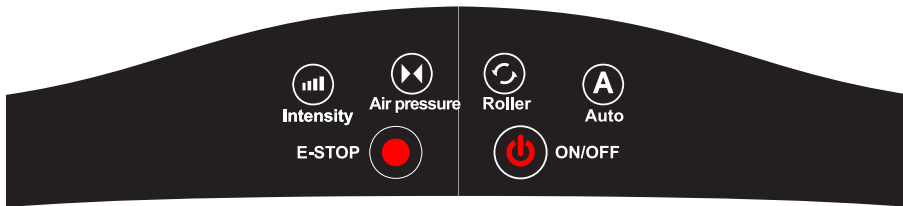Getting started:

Powering the machine , switch to the “-“ on the front box , underside of the right .

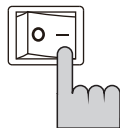Press the “  ” button , start into the program1 .


Descriptions :

3-auto programs are available for use , press this button will cycle to display .

 <b>Auto</b>	<p>Program1(red light) : Combining air pressure , roller , swing and acupressure massage functions .</p> <p>Program2(yellow light) : Distinctive swing feet massage , give all new experience .</p> <p>Program3(blue light) : Emphasizing on heels massage with air pressure and roller massage .</p>
--	---

※During using also can press the intensity button to adjust .

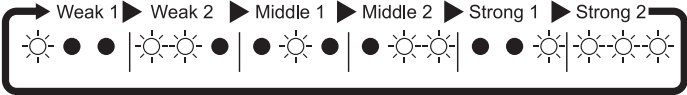 <p><b>Air pressure</b></p>	<p>Air pressure :</p> <p>※The initial setting is weak , press to adjust intensity .</p>  <p><b>Intensity</b></p>
 <p><b>Roller</b></p>	<p>Roller :</p> <p>Having a forward &amp; reversal rolling way .</p> <p>※The initial setting is slow , press to adjust speed.</p>  <p><b>Intensity</b></p>
 <p><b>Intensity</b></p>	<p>After activating the automatic system, the initial intensity is set at Weak 1. Press this button to adjust the intensity as following circulation lights:</p> <div style="text-align: center;">  <p>Blue    Yellow    Red</p> </div> 
<p><b>E-STOP</b></p> 	<p>Pause key :</p> <p>※All the motions would be stopped when press this button.</p>

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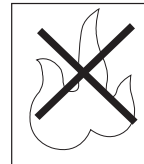
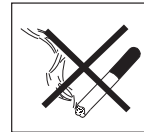
# CARE AND MAINTENANCE

## TROUBLESHOOTING

Please surveied the foot massager when you use it often

Please checked up the foot massager when you have a long time not in use.

SITUATIONS	SOLUTION
The foot massager doesn't work ,ON/OFF button don't light up	1.Check whether the plug is well connecting and working well 2.Check whether the switch on the foot massager is turnen to [I]positon
Any sound during work	Under those sound and feeling is normal phenomenon 1. The sound when rolling and air pressure massage 2. The different after adjust. 3. Internal sound when massage 4. Massage wheel grated on the cloth 5. Expansive sound and exhaust sound during air massage

## EMERGENCY TREATMENT

1. An exception occurred during the foot massager work,please turn off the switch and pull out the power line.
2. If the power line is broken,please notify our customer service to processing.

**Warning:** Except professional you don't dismantling,repair and change the foot massager ,it may reduce electric shock and risk of fire.

## SPECIFICATIONS

Model	FeelFoot (Series TF)
Rated Voltage and Frequency	<input type="checkbox"/> 100~120V <input type="checkbox"/> 220~240V <input type="checkbox"/> 50Hz
Rated Power	50W
N.W./G.W.	8.4/9.5kg
Dimension	L380xW450xH256(mm)
Data of manufacture:	look at the massager

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Авторское право © 2012г.
Ergonova GmbH